



Technology for Life Science

Non-mydriatic Retinal Camera

KOWA *nonmyd* 8

KOWA *nonmyd* 8s

**New**

Our experience for you







Technology for Life Science

# High resolution photography allows better clinical decisions!

## Featuring

- ▶ Non-mydratric color retinal photography [fig1]
- ▶ Fundus Autofluorescence (FAF) with low flash intensity\*\* [fig2]
- ▶ Anterior photography [fig3]
- ▶  $\phi$  3.3 mm minimum pupil diameter
- ▶ Multiple internal fixation targets
- ▶ Japanese manufactured optics for optimum performance

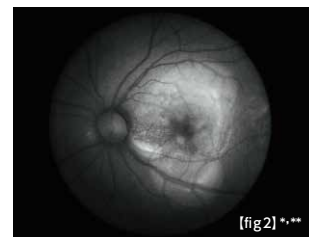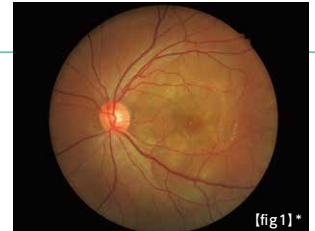

The KOWA nonmyd 8s is also available as basic model for color photography only.

\* Provided by Showa University  
\*\* Available by KOWA nonmyd 8

## Specifications

Photography modes	nonmyd 8 : Color, FAF nonmyd 8s : Color
Field angle	45°
Working distance	30 mm
Minimum pupil diameter	$\phi$ 4.0 mm (SP $\phi$ 3.3 mm)
Compensation range of examined eye	Without compensation : -12 D ~ +13 D Compensation + : +10 D ~ +35 D Compensation - : -32 D ~ -10 D
Focusing	Split luminous bars coincidence
Working distance adjustment	2 luminous dots indication type
Monitor	5.7 inch LCD

Light source	Observation : Near-infrared LED Photography : Xenon flash lamp
Internal fixation target	Central, Disc, Macula, mosaic 8 positions
Optical component adjustment range	Forward/Backward : 40 mm Leftward/Rightward : 98 mm Vertical (electric) : 27 mm
Chin rest adjustment range	55 mm (electric)
Power supply	Input : AC 100-240 V 50/60 Hz Power consumption : 150 VA
Dimensions	310 (W) × 504 (D) × 548 (H) mm
Weight	21 kg / 46 lbs (excluding the attached digital camera)

nonmyd is a registered trademark of Kowa Company, Ltd. in the United States, Germany, Japan and other countries.  
Specifications and appearances are subject to change without notice.

Distribution name : KOWA nonmyd 8  
Distribution name : KOWA nonmyd 8s



### Kowa Company, Ltd.

**World Sales Headquarters**  
4-14, Nihonbashi-honcho 3-chome, Chuo-ku,  
Tokyo 103-8433, Japan  
Phone: +81(3)3279-7639  
Facsimile: +81(3)3279-7541  
URL: <http://www.kowa.co.jp/e-life/>


**Kowa Company, Ltd. Chofu Factory**  
3-1, Chofugaoka 3-chome, Chofu,  
Tokyo 182-0021, Japan

### Kowa American Corporation

**Medical Devices Division**  
20001 South Vermont Avenue,  
Torrance, CA 90502, U.S.A.  
Phone: +1(310)327-1913 Facsimile: +1(310)327-4177  
URL: <http://kowa-usa.com/>

### Kowa Europe GmbH

Immermannstr. 43B, 40210 Düsseldorf, Germany  
Tel: +49 (0) 211-179354-18  
Fax: +49 (0) 211-161952  
Email: [medical@kowaeurope.de](mailto:medical@kowaeurope.de)  
Web site: [www.kowamedical.com](http://www.kowamedical.com)

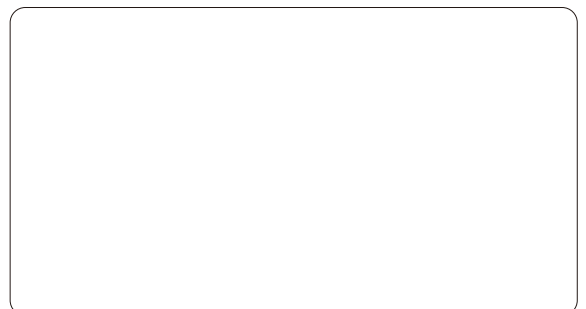

[REF] K9L61/63\_1699300CO  
Printed in Japan.

